

tecme

neumovent advance

Comprehensive monitoring  
Volumetric capnography  
High-flow oxygen therapy  
Spontaneous Breathing Test  
AVA (Adaptive Ventilatory Assist)

High-end mechanical ventilator  
that provides invasive and non-invasive  
ventilation, adaptable to any  
patient category.

**tecme**  
committed to life

## Comprehensive monitoring

It is possible to monitor the most important parameters, from the beginning of the therapy, with P/V Flex, as well as the intermediate phases towards weaning, with P0.1 and Pimax, among other tests and variables.

## Volumetric capnography

Integrated monitoring of exhaled CO<sub>2</sub> and display of volumetric curves over time, as well as measurements of physiological dead space.

## High-flow oxygen therapy

The addition of this method expands the use of the ventilator to the stages that are prior and after the mechanical ventilation itself, by providing non-invasive assistance to oxygenation.


## Spontaneous breathing test

This new tool allows a spontaneous breathing test to be performed before weaning, monitored in real time and completely safe with dedicated alarms.

## AVA (Adaptive Ventilatory Assist)

The AVA mode is an adaptive closed loop mode with optimal control that guarantees minute ventilation of the patient generating a ventilatory pattern with the minimal respiratory work.

- Ventilation for adult, pediatric and neonatal patients.
- Available for invasive and non-invasive ventilation.
- High-flow oxygen therapy.
- Automatic leak compensation.
- Integrated volumetric capnography software.
- Comprehensive respiratory mechanics module.
- Neonatal proximal flow sensor.
- Endotracheal or tracheotomy tube compensation.
- Intra-hospital transportation.
- Spontaneous Breathing Test (SBT).
- 72 hours of trend storage.
- Low cost maintenance.
- Built-in battery with capacity higher than 2.5 hrs.
- Numerical screen.
- 2-year warranty.

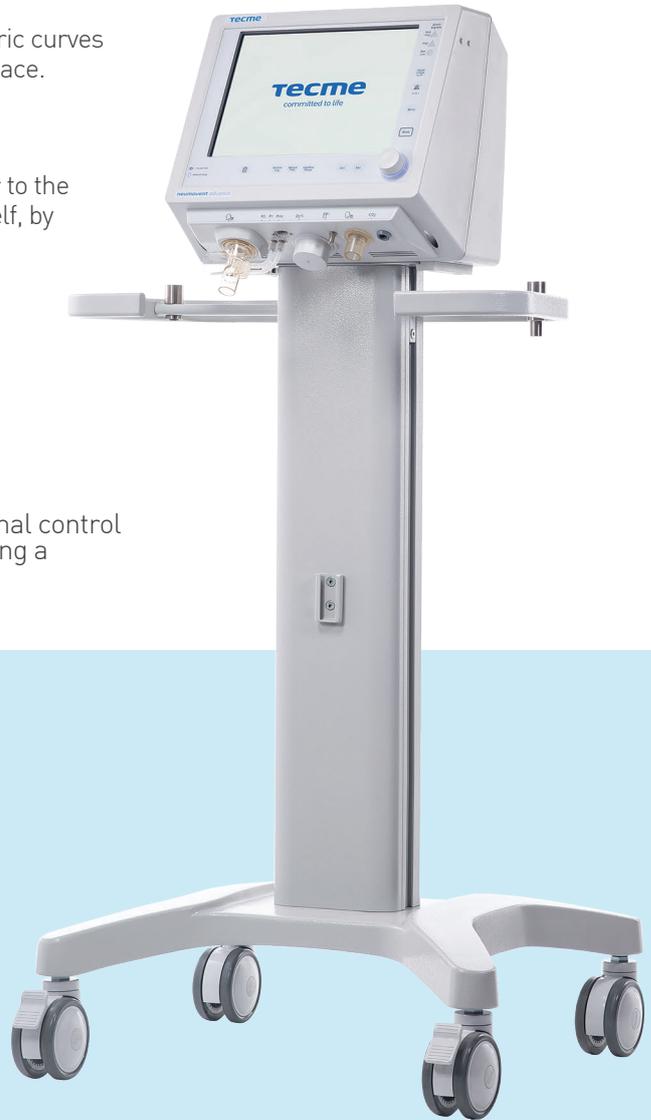

50 years of innovation and development in mechanical ventilators.

[www.tecmeglobal.com](http://www.tecmeglobal.com) / [info@tecmeglobal.com](mailto:info@tecmeglobal.com)

### Córdoba

Av. Circunvalación (N-O) - Agustín Tosco 3040  
X5008HJY - Córdoba, Argentina

P: +54 351 4144600 - F: +54 351 4144605

### Buenos Aires

Trole 557 - 2º Piso, Nueva Pompeya, C.A.B.A.,  
C1437DKK - Buenos Aires, Argentina

P: +54 11 1561387326

### Atlanta

2825 Pacific Drive, Suite B  
Norcross, GA 30071 - USA

P: +1 770 409 9172 - F: +1 770 729 8176